



## The FORMACS<sup>LT</sup> TOC Analyzer



The measurement of total organic carbon (TOC) constitutes an essential part of many water quality studies. Using a combination of UV light and wet oxidation the Formacs<sup>LT</sup> TOC Analyzer is particularly suitable for the analysis of complex sample matrices in a wide variety of applications including; drinking, ground, sea, surface and waste waters, effluents, soil digests, boiler feed water, etc.


The principle components of the Formacs<sup>LT</sup> TOC Analyzer include an internal sparging system, continuous UV-promoted persulphate digestion, on-line stripping and NDIR (Non Dispersive Infrared) detection. The latest version Windows<sup>®</sup> based software fully compliments this set up for all data handling and presentation.

The operator can choose to include or exclude inorganic and volatile compounds. The process of removing these is automatically performed by acidifying and sparging the sample with nitrogen gas. With addition of potassium persulphate and exposure to UV light, the organic compounds of the sample are oxidized to CO<sub>2</sub>. The CO<sub>2</sub> is expelled from the solution by a stripper unit and is measured quantitatively using NDIR detection.

The instrument includes a fully integrated Random Access Sampler. The sampler has a capacity of 80 sample positions. A contamination free top stirrer is fitted as a standard feature.

### Features:

- Measures TC, TIC, TOC, DOC and NPOC
- Highly sensitive (ppb levels)
- Wide range of sample matrices can be handled
- Minimum bench space required
- Fully computer controlled
- Storage of multi range calibration curves
- Real time graphics
- Easy to use software
- Methods according to Standard Methods, EPA, CEN, DIN, ISO etc.

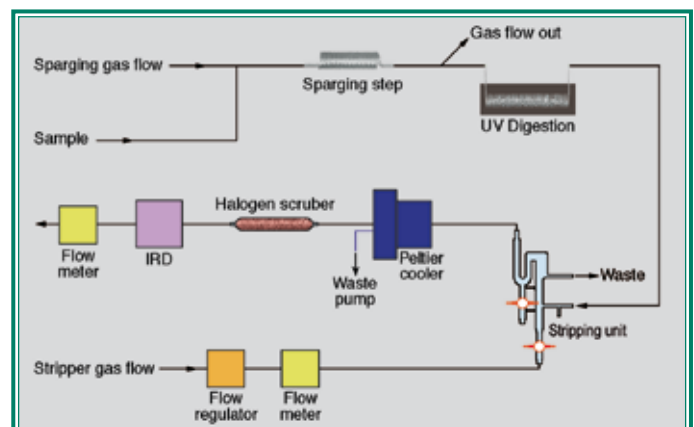

Flow diagram of the Formacs<sup>LT</sup>

**Sample****Oxidation efficiency**

Glycine	97.9 %
EDTA	94.9 %
Glucose	100.0 %
Caffeine	104.4 %
4-Aminoantipyrine	99.8 %
Thiamine	97.3 %
Barbituric acid	100.8 %
Sucrose	101.1 %
Potassium hydrogenphthalate	100.0 %

**General Characteristics**

Analytes	Total Carbon (TC), Total Inorganic Carbon (TIC), Total Organic Carbon (TOC), Dissolved Organic Carbon (DOC), Non Purgeable Organic Carbon (NPOC)
Method	UV promoted persulphate digestion with synchronous dual wavelength Infra Red detection (NDIR)
Samples	Drinking-, ground-, waste-, sea-, cooling-, surface-, pharmaceutical waters
Sample introduction	Injection by peristaltic pump
Features	Integrated random access 80 position sampler Extendable with Primacs <sup>MCS</sup> – Add-on module for TOC in solid samples Contamination free top-stirring Complies with international regulations such as EPA 415.2, DIN 38409 H3 Standard Methods 5310C

**Operational and Performance Characteristics**

Measuring range	Up to 2 ppm C - Up to 10 ppm C - Up to 100 ppm C - Up to 1000 ppm C
Detection limit	0.02 ppm C (lower ranges can be achieved but depend on laboratory conditions)
Reproducibility	Ranges < 2 ppm within 3% of f.s. range. - Ranges > 2 ppm within 2% of f.s. range
Analysis Sequence	Continuous, results are displayed every 4-5 minutes
Data processing	Detector signal: Area calculation. Data: Multi point linear regression, exclusion of results and recalculations, statistical calculation

**Physical Characteristics**

Gas Combustion	Carrier gas: Nitrogen or purified compressed air.
Power requirements	110 – 120 V, 220 – 240 V, 60/50 Hz
Dimensions (hxdxw)	58 x 53 x 72 cm
Weight	53 kg

**Skalar's Headquarters**  
**Skalar**  
 P.O. Box 3237  
 4800 DE Breda  
 The Netherlands  
 Tel. +31 (0)76 5486 486  
 Fax. +31 (0)76 5486 400

**USA**  
**Skalar, Inc.**  
 5995 Financial Drive, Suite 180  
 Norcross, GA 30071  
 Tel. + 1 770 416 6717  
 Toll Free: 1 800 782 4994  
 Fax. + 1 770 416 6718

**United Kingdom**  
**Skalar (UK) Ltd.**  
 Breda House,  
 Millfield Industrial Estate,  
 Wheldrake, York, YO19 6NA  
 Tel. + 44 (0)1904 444800  
 Fax. + 44 (0)1904 444820

**Germany**  
**Skalar Analytic GmbH**  
 Gewerbestraße Süd 63  
 41812 Erkelenz  
 Germany  
 Tel. + 49 (0)2431 96190  
 Fax. + 49 (0)2431 961970

**Austria**  
**Skalar Analytic GmbH**  
 Am Anger 22  
 A-7451 Oberloisdorf  
 Austria  
 Tel. + 43 (0)2611 2023411  
 Fax. + 43 (0)2611 2023412

**France**  
**Skalar Analytique S.A.R.L.**  
 79, Avenue Aristide Briand  
 94110 Arcueil  
 France  
 Tel. + 33 (0)1 4665 9700  
 Fax. + 33 (0)1 4665 9506

**Belgium**  
**Skalar Belgium bvba**  
 Antwerpsestraat 126  
 2850 Boom  
 Belgium  
 Tel. + 32 (0)3888 9672  
 Fax. + 32 (0)3844 3441



Email: [info@skalar.com](mailto:info@skalar.com)  
 ©Copyright Skalar 2008

Internet: [www.skalar.com](http://www.skalar.com)  
 Publication No. 0406003F. US

Skalar reserves the right to change specifications and appearance of the equipment without further notification.